



Fruit Texture Analyser

Get firm - insist on doing it electronically

The Fruit Texture Analyser (FTA) measures fruit firmness at the push of a button and captures the results to a PC. Tests are conducted at standard depths and speeds of penetration, ensuring **accurate** and **repeatable** results.

The FTA is ranked best fruit firmness testing instrument by a [Washington State University](#) study. ([Download full report](#)).

Get a taste for accurate measurement . . . every time!

Test apples, pears, peaches, nectarines and kiwi-fruit - with **quick** , **easy** **reliable** results.



Some firm advantages

- The FTA connects to a standard PC and the software runs on Windows 98/NT4/ME/2000/XP.
- The test results can be processed as follow:
 - Save the results to disk (or network) in ASCII format.
 - Capture the test results with your own software directly from the FTA using the control functions (supplied with the system in a DLL).
- The FTA can directly be accessed from Excel, Word or Access. Click [here](#) to see how it works from Excel. Samples are included on CD.
- Installation only takes a few minutes.
- The FTA requires minimum maintenance - only cleaning to remove juice.
- The FTA can also be controlled from:
 1. [ART-Soft database](#) which forms part of the ART System (Apples Ripening Test). The software can be purchased directly from [UP GmbH](#), Germany.
 2. FQMS database (Fruit Quality Management System). Contact GÜSS or [Hortec](#) for further information on this software.

Also measure size and weight


The EFM (Electronic Fruit Size Measure) can be connected to the FTA to measure fruit diameter electronically. Each measured size will be immediately available on the computer. Fruit diameter from 12 mm to 100 mm can be measured.

Note: the EFM is an optional attachment.



A Scale attachment can also be connected to the FTA to measure fruit weight electronically. The Scale is available in capacities ranging from 3 kg to 35 kg.

Note: the Scale is an optional attachment.

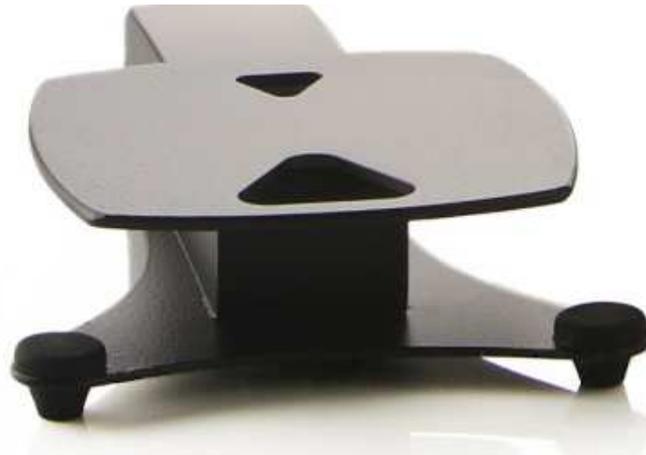

LCD/Keypad

With the LCD/Keypad attached the FTA can be operated without a connection to a PC. All the relevant measurements - firmness, size, minimum, maximum, average etc., will be displayed on the LCD.

Note: the LCD/Keypad is an optional attachment



Features and benefits


- The FTA can be controlled with **your own software** . The control functions can be accessed through a DLL file. Integrating the FTA with existing QA software systems should not take longer than two hours of programming. [Download SDK](#)
- The electronic measure system ensures **accurate, reliable firmness readings** .
- The motorized system **eliminates operator discrepancies** and ensures accurate test distance (depth).
- Electronic speed control ensures **accurate probe speed** and also ensures that tests are carried out under standard conditions.
- Built-in safety features protect against overload.
- Software can be configured to allow for customised output data and enables the FTA data to be incorporated into existing quality control systems.
- All operating parameters can be regulated via the control PC (probe speed, calibration, threshold load, depth).
- **Low maintenance** - NO lubrication or adjustments are required.
- Software to generate reports from the recorded data is included in the package .
- Measuring units:


Firmness: kg, pounds, Newton.
Depth: millimeter, inch.

It's all in the system

To run the control software, you need a computer with at least:

- Windows XP,7 or 8.1
- USB or Serial Port (DB9 connector)
- CD-ROM drive

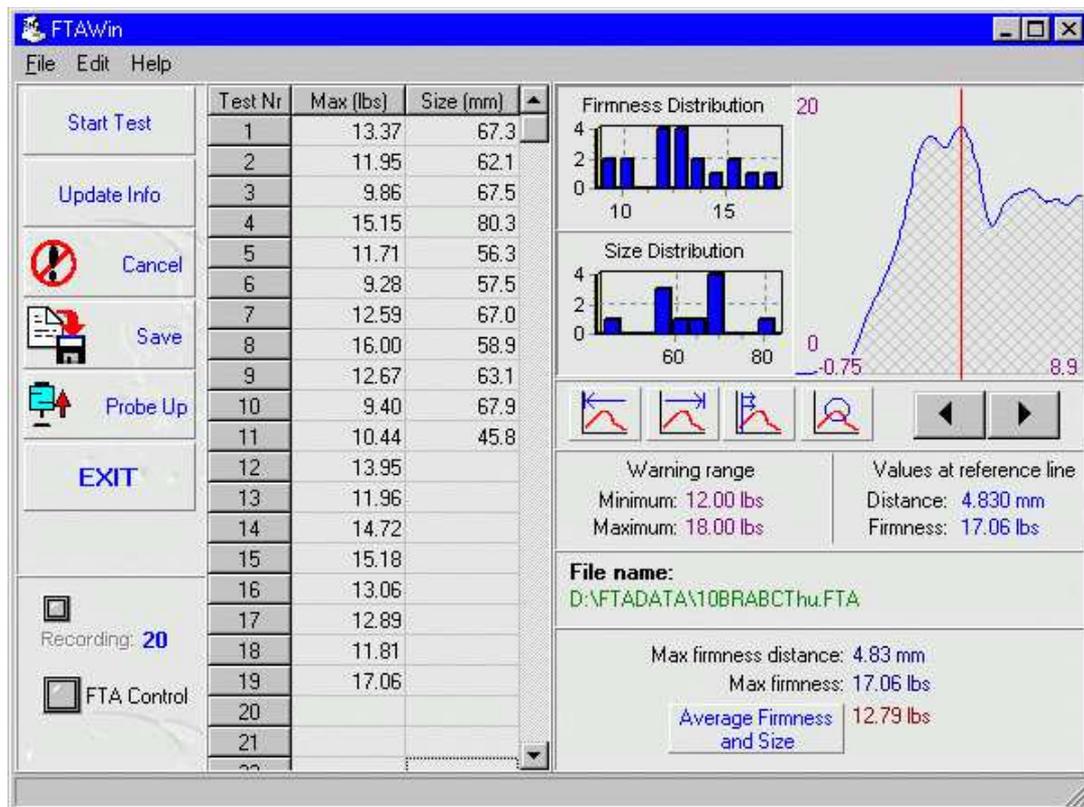

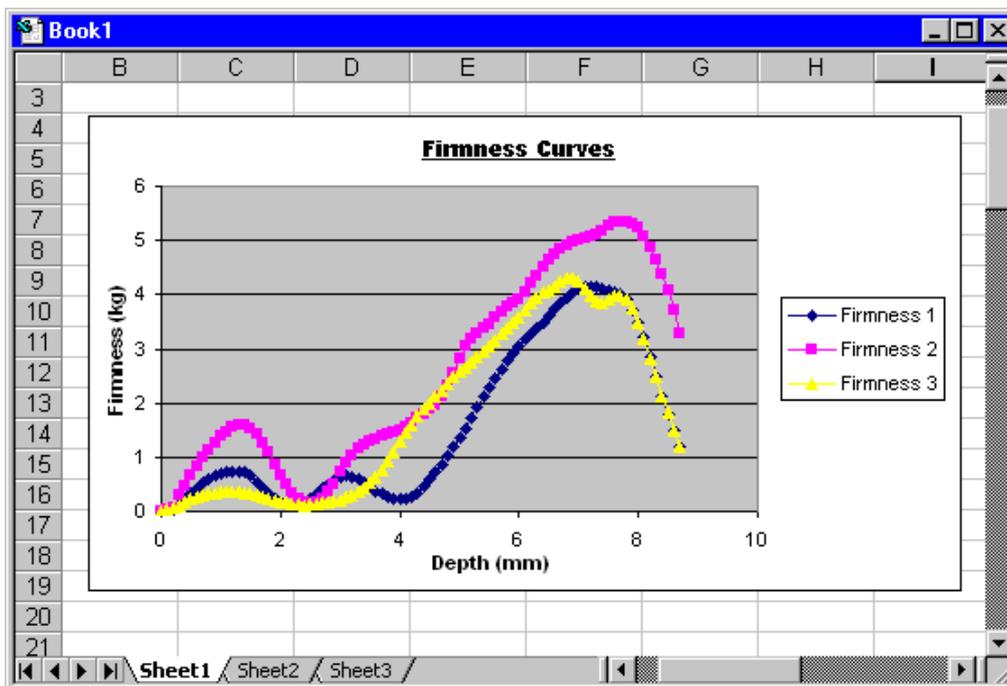

The firmness data can easily be copied to Excel to create graphs for reports or presentation.

TECHNICAL DATA

	GS-15	GS-25
Speed of probe	1 ~ 40 mm/sec	1 ~ 40 mm/sec
Operating range	50 gram ~ 15 000 gram	50 gram ~ 25 000 gram
Maximum fruit size	120 mm	120 mm
Mass	9 kg (approx.)	9 kg (approx.)
Power consumption	50 watt (1)	50 watt (1)
EFM (size measure) port (2)	standard	standard
Scale attachment port (3)	standard	standard
LCD/Keypad port (4)	standard	standard
Digital caliper port (5)	standard	standard
Supply voltage	30 V DC	30 V DC
USB	Yes	Yes

- 1) under max. load.
- 2) does not include EFM attachment
- 3) does not include scale attachment
- 4) does not include LCD/keypad attachment
- 5) connects to *Mitutoyo* calipers only

[[Home](#)] [[Top](#)]

home page www.gusstoday.com 